



# Combivent Respimat

Combivent Respimat® (albuterol & ipratropium bromide)

100mcg/20mcg



Medicine	Need a	My dose
Quick Reliever (soft mist inhaler)	No	Check Asthma Action Plan

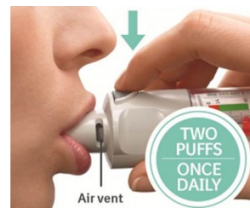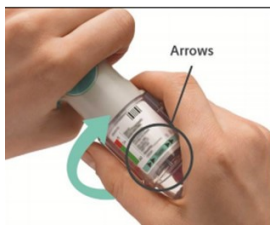

1. Use the TOP Acronym
2. T- TURN the clear base in the direction of the arrows on the label until it clicks (keep the cap closed while doing this)
3. O- OPEN the cap until it snaps fully open
4. P- Breathe out slowly away from device, close your lips around the mouthpiece without covering air vents. While taking a slow deep breath through your mouth, PRESS the dose release button and continue to breathe in slowly. Hold breath as long as you are comfortable.
5. Repeat steps for each puff
6. Close the cap until you use your inhaler again



Use the camera on your phone to scan this QR Code for the Combivent Respimat Website:

1/6/2021

There is information about:

- ◆ How to use the Respimat
  - ◆ Click on “Preparing your Respimat inhaler” and “Using Respimat every day” to watch videos